

The **PRECISION** you
require with the
VERSATILITY you need
to be successful



Radiofrequency Energy Sources & Accessory Products




Dermatology
Ear, Nose, and Throat (ENT)
Facial Plastic
General Surgery
Gynecology
Ophthalmology / Oculoplastics
Oral / Maxillofacial Surgery
Plastic Reconstructive
Primary Care
Podiatry
Proctology

ellman[®]

Experts in Precision Surgery

Established in 1959

Radiosurgical Energy Sources	4
Five Distinct Waveforms	4
Surgitron® 4.0 MHz Dual RF™	5
Surgitron® Radiolase® I	6
Surgitron® Radiolase® II	7
Surgitron® FFPF EMC	8
Radiosurgical Accessories	9
Neutral Plates	9
Footswitches	10
Cables	11
Surg-e-Doc™	12
EMC Device Accessories	13
Handpieces	15
Bipolar Forceps	17
Bayonet Bipolar Forceps	19
ACe-Tip™ Non-Stick Forceps	20
RF SurgiPens™	22
RF Electrodes	23
Sterile/Single-Use Electrodes	24
Empire® microIncision™ Needles	25
Reusable Electrodes	27
ACe-Tip™ Electrodes	32
Reusable Electrode/ACe-Tip™ Sets	33
Electrode Accessories	35
Educational Materials	36
Compact Discs	36
Patient Education Brochures	36
Textbooks	36
Smoke Evacuation System	37
Accessories	38
Multi Specialty Accessories	39
Otorhinolaryngology (ENT)	40
General Surgery/Proctology	47
Gynecology	48
Plastic Surgery	50
Dermatology	52
Non Ablative Skin Tightening Handpieces	53
Ophthalmology/Oculoplastics	54
Family Practice/Podiatry	56
Hemostasis Electrode	57
Index	58

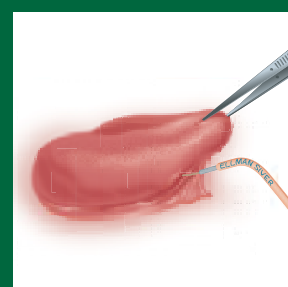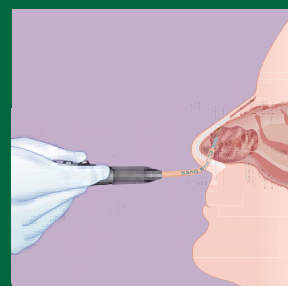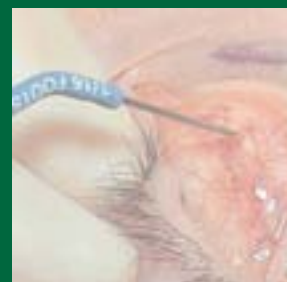


RADIOSURGICAL ENERGY SOURCES

5 Distinct Waveforms

1


Fully Filtered (Cut)

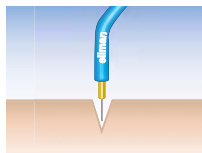

- Micro-smooth cutting
- Negligible lateral heat
- Minimal cellular destruction
- Best cosmetic results. Fastest healing^{1,2}
- Ideal for skin incision and biopsy
- 4.0 MHz

2



Fully Rectified (Cut/Coag)

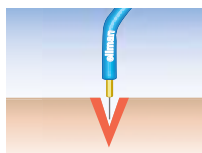

- Cutting with hemostasis
- Ideal for sub-cutaneous tissue dissection and planing. Especially useful in vascular areas while producing minimal amounts of lateral heat and tissue damage
- 4.0 MHz

3


Partially Rectified (Hemo)

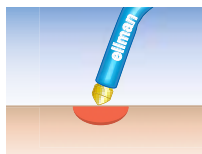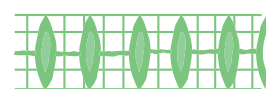

- Coagulation / Shrinkage
- Hemostasis with controlled penetration
- Ideal for cutting with hemostatic control
- 4.0 MHz

4

Fulguration

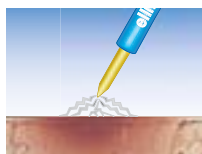


- Maximum penetration and hemostasis
- Ideal for intentional tissue destruction
- 4.0 MHz

5

Bipolar

- Pinpoint, micro-coagulation
- Minimal tissue adherence to forceps
- Minimal charring or tissue necrosis
- Ideal for coagulation in and around critical anatomy
- 1.7 MHz

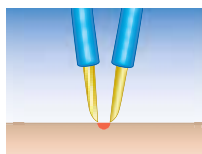


Clinical Citations

1. Bridenstine, J.B., Derm Surgery (1998); vol 24, p397-400
2. Aferzon, M. Derm Surgery (2002); vol 28, p735-738

ellman's radiofrequency technology has received worldwide acclaim and recognition from leading surgeons, clinics, universities and authors.

TRANSMISSION ELECTRON MICROSCOPY SHOWING LATERAL THERMAL DAMAGE ON HUMAN OVIDUCTS BY DIFFERENT INSTRUMENTS

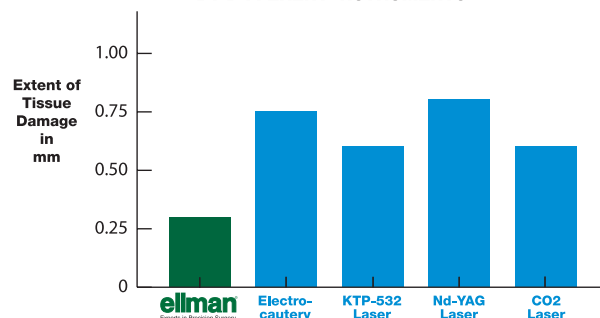

Reference - Oilvar, AC, et al, Ann Clin Lab Sci. 1999 Oct-Dec; 29(4): p281-5.

The ellman® radiowave energy source represents advanced radiowave technology that provides outstanding surgical control, precision and versatility. Clinical benefits include reduced post-operative discomfort, minimal scar tissue formation, enhanced healing and excellent cosmetic results.

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

RADIOSURGICAL ENERGY SOURCES

Surgitron® Dual RF™



Surgitron® Dual RF™/120 IEC High Frequency — Low Temperature RF Energy Source

The patented Surgitron® Dual RF™ Frequency unit represents advanced radiowave technology that provides unparalleled surgical control, precision and versatility. The high frequency of 4.0 MHz minimizes heat dissipation and thus cellular alteration. The system generates 120 watts of power and operates at 4.0 MHz in a monopolar mode and 1.7 MHz in bipolar mode.

Clinical Benefits Include:

- Reduced post-operative discomfort¹
- Minimal scar tissue formation²
- Maximum readability of histologic specimen³
- Enhanced healing⁴
- Excellent cosmetic results⁵

Clinical Citations

1. Ericsson, E et al, The Laryngoscope (2007); vol 117, p654
2. Botero, G.E.S, J Otol Head Neck Surgery (1996); vol 24 (1), p69
3. Silverman, E.B et al, Veterinary Surgery (2007); vol 36, p50-56
4. Bridenstine, J.B., Derm Surgery (1998); vol 24, p397-400
5. Aferzon, M, Derm Surgery (2002); vol 28, p735-738

Patented Surgitron® 4.0 MHz
Dual RF™/120 IEC

***Item # IEC3A-S30**

Shipped complete with:

DF-FSC#	Dual Footswitch & Cable
IEC3A-SL1	Instruction Manual
IEC-HP1	IEC Foot-Controlled Handpiece
IEC-3FHPB	IEC 3-Button Fingerswitch Handpiece
H79	Handpiece Clip
VCD/PC5IEC3	Instructional CD ROM
IEC-BC/I	3- Bipolar Sterile/Single-Use Cables
IEC-NPC	IEC Neutral Plate
IEC-3FHPB/D1	Sterile/Single-Use 3-Button Fingerswitch Handpiece
*IEC-PC110 (for 110v)	power cord
*IEC-PC220 (for 220v)	power cord
DSEP40	40 Medical Electrodes Assortment Package Sterile/Single-Use
"J" Series	Bipolar Forceps (Item number determined by selection) see page 14 & 15

**(Specify voltage/country specific
plug available
110-220V, EU, UA, ANZ)*

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

RADIOSURGICAL ENERGY SOURCES

Surgitron® Radiolase® I



Surgitron® Radiolase® I High Frequency — Low Temperature RF Energy Source

The **ellman®** Surgitron® Radiolase® I cuts and coagulates soft tissue using monopolar energy. This system generates 50 watts of power and operates at 4.0 MHz.

Patented Surgitron® Radiolase® I

***Item # RL-S50A**

Shipped complete with:

IEC-HP1	Foot-Controlled Handpiece
IEC-NPC	Neutral Plate
RL-SL1	Operation Manual
VCD/PC5RRL	Video
RL-FSC	Footswitch
H79	Clip
IEC-1FHPA	One-Button Fingerswitch Handpiece
*IEC-PC110	(for 110v) Power Cord
*IEC-PC220	(for 220v) Power Cord
DSEP40	40Medical Electrodes Assortment Package Single-Use

**(Specify voltage/country specific
plug available
110-220V, EU, UA, ANZ)*

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

RADIOSURGICAL ENERGY SOURCES

Surgitron® Radiolase® II



Surgitron® Radiolase® II High Frequency — Low Temperature RF Energy Source

The **ellman**® Radiolase® II affords physicians the ability to utilize both monopolar and bipolar modalities.

This system generates 50 watts of power and operates at 4.0 MHz in both monopolar and bipolar modalities.

Patented Surgitron®
Radiolase® II

***Item # RL2-S50**

Shipped complete with:

IEC-HP1	Foot-Controlled Handpiece with Locking Chuck
IEC-NPC	Neutral Plate
RL2-L1	Operation Manual
RL2-FSC#	Dual Footswitch
IEC-1FHFA	One-Button Fingerswitch Handpiece
IEC-BC/I	Bipolar Cable
*IEC-PC110	(for 110v) Power Cord
*IEC-PC220	(for 220v) Power Cord
DSEP40	40 Medical Electrode Assortment Package Sterile/Single-Use
"J" Series	Bipolar Forceps (Item number determined by selection) see page 14 & 15

**(Specify voltage/country specific
plug available
110-220V, EU, UA, ANZ)*

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

RADIOSURGICAL ENERGY SOURCES

Surgitron® FFPF EMC



Surgitron® FFPF EMC High Frequency — Low Temperature RF Energy Source

The **ellman**® Surgitron® FFPF EMC cuts and coagulates soft tissue using monopolar energy. This system generates 90 watts of power and operates at 3.8 MHz

Patented Surgitron® FFPF EMC Device

***Item # EMC 110A** (for 110v) Power Cord

***Item # EMC 220A** (for 220v) Power Cord

Shipped complete with:

HP1	Foot-Controlled Handpiece with Locking Chuck
H4	EMC Neutral Plate
*PC	EMC Power Cord
DSEP40	40 Medical Electrode Assortment Package Sterile/Single-Use
EMCL1	Operating Manual
H79	Handpiece Clip
DVD/EMC1	Instructional Video

**(Specify voltage/country specific
plug available
110-220V, EU, UA, ANZ)*

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

RADIOSURGICAL ACCESSORIES

Neutral Plates

IEC Neutral Plate

IEC-NPC


Reusable, coated
For Dual RF™ and Radiolase® Devices
(1 pc)



IEC Disposable Neutral Plate-Standard Size

IEC-NPD

Single-Use
For Dual RF™ and Radiolase® Devices
(25 pcs)

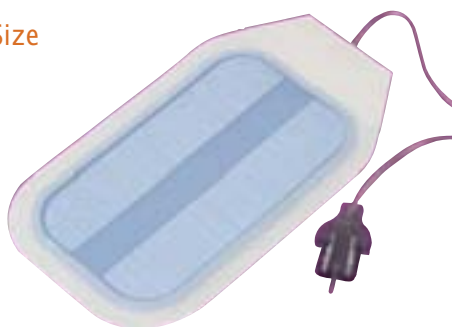

IEC Disposable Neutral Plate-Small Size

IEC-NPD-S

Single-Use
For Dual RF™ and Radiolase® Devices
(25 pcs)



EMC Neutral Plate

H4

Reusable, coated
For the EMC Device
(1 pc)



With the **ellman®**
Reusable Neutral Plate,
no skin contact required.

The **ellman®** Single-Use Neutral Plate is a sensing, pre-gelled, self-adhering, dispersive pad. The effective electrodes are two aluminum oil conductors, covered by a hydrogel conductive adhesive. The conductive hydrogel provides a low-current density return path from the patient to the RF Device while the two foils provide a means of sensing contact quality.

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

RADIOSURGICAL ACCESSORIES

Footswitches

Dual Frequency Footswitch and Cable

DF-FSC#

For Surgitron® Dual RF™ Device
(1 pc)



Single-Action Bipolar Footswitch and Cable

DF-FSC/3

For Surgitron® Dual RF™ Device
(1 pc)



Radiolase I Footswitch and Cable

RL-FSC

For Surgitron® Radiolase® I Device
(1 pc)



Radiolase II Monopolar/Bipolar Footswitch

RL2-FSC#

For Surgitron® Radiolase® II Device
(1 pc)



Pneumatic Footswitch

SVFSC

For use with Surg-e-Vac™
Smoke Evacuation System
(1 pc)



For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

RADIOSURGICAL ACCESSORIES

Cables

Bipolar Cable

IEC-JX/B

Reusable

For use with Surgitron® Dual RF™ or Radiolase® II Devices (1 pc)

IEC-BC

Sterile/Single-Use

Minor Surgery Bipolar cable

For use with Surgitron® Dual RF™ or Radiolase® Devices (5 pcs)

JX/E

Reusable

For use with EMC Device (1 pc)

Quick-Connect™ Cable

For use with SurgiPen™ Suction Tonsil and Suction Adenoid Probes, Disposable Suction Coagulator, Mono Stan RF Forceps, Rodin™ Micro Surgery handpiece, Micro Fibre® Electrode.

IEC-SP020

Reusable

For use with Surgitron® Dual RF™ or Radiolase® II Device (1 pc)

SP020

Reusable

For use with EMC Device (1 pc)

Snare Cable

For use with Polyp Snare™

IEC-H7B

Reusable

For use with Surgitron® Dual RF™ or Radiolase® II Device (1 pc)

H7B

Reusable

For use with EMC Device (1 pc)


Lap/Endo Monopolar Cable

IEC-UC

Reusable

For use with Surgitron® Dual RF™ or Radiolase® II Device (1 pc)

Monopolar cable for Endoscopic instruments

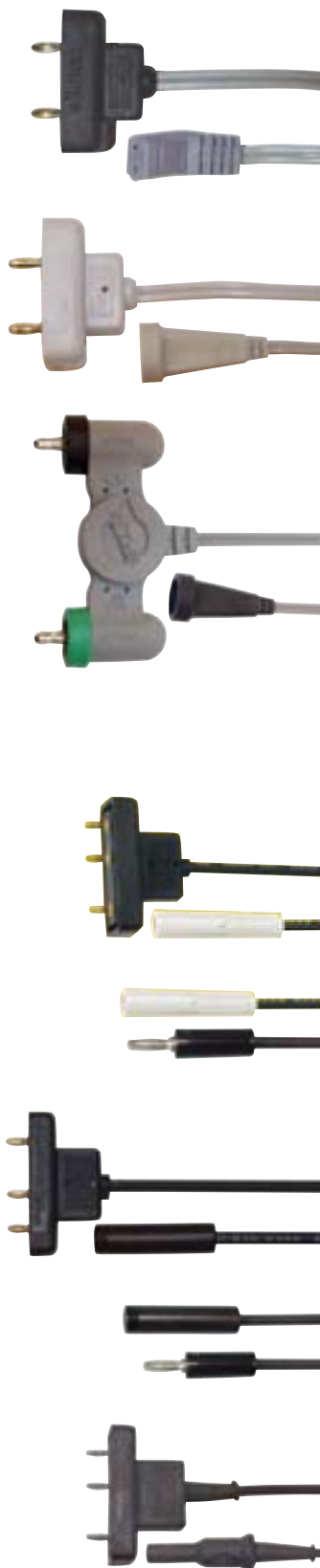

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

RADIOSURGICAL ACCESSORIES

Surg-e-Doc™



*Surg-e-Doc™ shown with
Surgitron® Dual RF™*

Surg-e-Doc™ SDS30

The Surg-e-Doc™ offers the following advantages:

Extra Monopolar port

The Surg-e-Doc™ offers the ability to utilize an extra monopolar port. This added versatility eliminates the need to change monopolar electrodes within the same handpiece during a procedure.

Simultaneous activation of Surg-e-Vac™

The Surg-e-Doc™ extends the life of the vacuum unit by automatically activating the Surg-e-Vac™ only when the RF energy source is activated. Moreover, the Surg-e-Vac™ will continue vacuuming for 5 seconds after the RF energy source is deactivated. This allows for automatic removal of residual plume.

Less Clutter

The Surg-e-Doc™ eliminates the need for separate footswitch activation of Surg-e-Vac™. This reduces clutter since the only footswitch needed is that of the RF energy source.

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

RADIOSURGICAL ACCESSORIES

Adapters for Surgitron® EMC

Universal Fingerswitch Adapter

EMCKA

Reusable

For use with EMC Device, handpiece not included
(1 pc)



Patented 1-Button Fingerswitch Handpiece with Locking Chuck

XK01

Reusable

For use with EMC Universal Fingerswitch Adapter
(2 pcs)



EMC Universal
Fingerswitch Plug

The **ellman®** Universal Fingerswitch Adapter plugs into the side of your Surgitron® FFPF EMC. The adaptor eliminates the need for a foot pedal. The standard foot-controlled handpiece is still operational.

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

RADIOSURGICAL ACCESSORIES

For Surgitron® EMC

EMC Steri Caps™

EMCH35

Reusable
Steri Cap™ (3 for large dial and
3 for small dial) for use with
EMC only, 3-Blue and 3-Black
(6 pcs)



Steri Caps™ enable you to
control in total sterility with
autoclavable protective
dial covers.

Patented Bipolar Converter

EMCH95A

For use with EMC Device
(1 pc)



The Bipolar Converter
provides the surgeon the
convenience of switching the
ellman® Surgitron® EMC from
the Handpiece in monopolar
mode to the Bipolar Forceps
mode with the simple flick of a
switch. This eliminates the
need to plug and unplug the
cords during surgery.

Bipolar Converter Steri Caps™

H96

Caps for activation switch
for use with Bipolar Converter
(3 pcs)



Wall Bracket

HH

Wall Bracket with hardware.
For use with EMC only
Swivel arm and
swivel plate.
(1 pc)



The Wall Bracket space
saving accessory attaches
easily to any vertical surface.
The **ellman® Wall Bracket** will
support up to 40 lbs. and is
ideal for the Surgitron® EMC.

For orders and/or pricing information, please call **ellman®** at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

RADIOSURGICAL ACCESSORIES

Handpieces



Patented 3-Button Fingerswitch Handpiece with Locking Chuck

IEC-X3FHPB

Reusable (2 pcs)

Accepts 1/16" and 3/32" Shaft Electrodes

For Surgitron® Dual RF™



Patented Sterile 3-Button Fingerswitch Handpiece with Locking Chuck

IEC-3FHPB/D

Sterile/Single-Use (10 pcs)

Accepts 1/16" and 3/32" Shaft Electrodes

For Surgitron® Dual RF™




Patented 3-Button Fingerswitch Blade Handpiece with Locking Chuck

IEC-XBHP3

Reusable (2 pcs)

Accepts Sterile/Single-Use Blades (assorted surgical blades including #10, and #15) and 1/16" Electrode

For Surgitron® Dual RF™

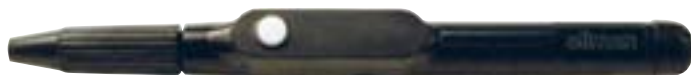

Patented 1-Button Fingerswitch Handpiece with Locking Chuck

XK01

Reusable (2 pcs)

Accepts 1/16" and 3/32" Shaft Electrodes

For Surgitron® EMC Device

The **ellman®** handpiece is ergonomically designed to fit comfortably in your hand therefore maximizing productivity by reducing operator fatigue and discomfort.

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

RADIOSURGICAL ACCESSORIES

Handpieces




Patented 1 Button Fingerswitch Handpiece

IEC-1FHPA

Reusable (2 pcs)

Accepts 1/16" and 3/32" shaft electrodes

One button fingerswitch handpiece for Radiolase® Devices

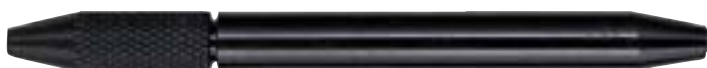

microInsulated Needle Handpiece with Locking Chuck

IEC-XH136 — For Surgitron® 4.0 Dual RF™ and Radiolase® II Devices (2 pcs)

XH136 — For Surgitron® EMC Device (2 pcs)

Reusable (1 pc). Smaller orifice for insertion of **ellman®** microInsulated needles only

Included in Cosmetic Treatment package



Blade Handpiece with Locking Chuck

IEC-H15 — Accepts 1/16" Shaft. For Surgitron® Dual RF™ and Radiolase® Devices

H15 — Accepts 1/16" Shaft. For Surgitron® EMC Device

IEC-QH15 — Accepts 3/32" Shaft. For Surgitron® Dual RF

QH15 — Accepts 3/32" Shaft. For Surgitron® EMC

Reusable (1 pcs). Foot Pedal activation

For use with assorted surgical blades including #10, and #15



Foot-Controlled Handpiece with Locking Chuck

IEC-XHP1 — For Surgitron® Dual RF™ and Radiolase® Devices

XHP1 — For Surgitron® EMC Device

Reusable (2 pcs). Accepts 1/16" and 3/32" shaft



Rodin® Microsurgery Handpiece

XRHP — Length (3.895") 98.93mm, Width (.460") 11.68mm

Reusable (2 pcs). To be used with IEC Quick-Connect™ Cable

Accepts 1/16" Shaft. For Surgitron® Dual RF™ and Radiolase® Devices


(See page 11 for Quick-Connect™ Cable, Cat. No. IEC-SP020)

(See page 33 for Rodin™ Microsurgery ACe-Tip™ Electrode Set, Cat. No. RACE7)

RODIN® 1

Rodin® RF Microsurgery Set

One RHP Rodin® Microsurgery Handpiece, one IEC-SP020 Quick-Connect™ Cable, one RACE7 Rodin® Microsurgery ACe-Tip™ Set, and instructions for use.

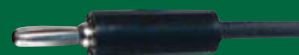

Surgitron® EMC
(Monopolar)



Surgitron® Dual RF™ and
Radiolase® (Monopolar)

Separate plugs are used in accordance with the corresponding device, the Surgitron® Dual RF™, Radiolase® and the Surgitron® EMC.

For orders and/or pricing information, please call **ellman®** at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

RADIOURGICAL ACCESSORIES

Bipolar Forceps



J1 Length 9.5cm, 0.4mm, Straight Tip



J1B Length 9.5cm, 0.4mm, Bent Tip



J2 Length 9.5cm, 0.3mm, Curved Flat/Angled Tip



J3 Length 7.8cm, 0.4mm, Angled Tip



J4 Length 7.8cm, 0.3mm, Straight Tip




J5 Length 9.2cm, 0.3mm, Bent Tip



J8 Length 16.3cm, 1.0mm, Straight Tip

ellman® Bipolar Forceps provide clean wet-field coagulation.

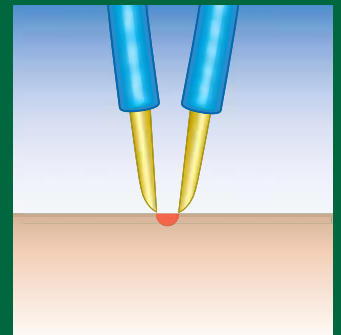

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

RADIOSURGICAL ACCESSORIES

Bipolar Forceps



J9 Length 16.2cm, 1.0mm, Straight Tip



J10 Length 16.2cm, 2.0mm, Straight Tip



J11 Length 12.2cm, 0.4mm, Straight Tip




J12 Length 10.9cm, 1.0mm, Straight Tip



J14 Length 10.8cm, 0.7mm, Straight Tip



J19 Length 13.6cm, 2.0mm, Straight Tip

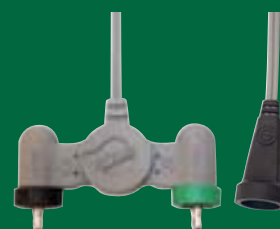

EMC Bipolar Cord
(see page 8)




Surgitron® Dual RF™ and
Radiolase® Bipolar Cord
(see page 8)

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

RADIOSURGICAL ACCESSORIES

Bayonet Bipolar Forceps

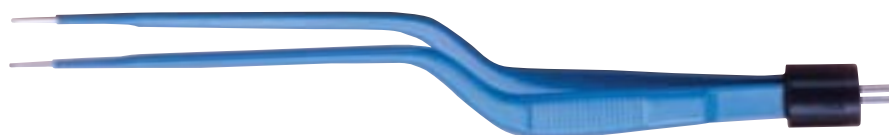


J7 Length 17.3cm, Tip Wide 1.0mm, Bayonet, Flat Tip

J7U Length 17.3cm, Tip Wide 1.0mm, Bayonet, 45° Up, Flat Tip


J7D Length 17.3cm, Tip Wide 1.0mm, Bayonet, 45° Down Flat Tip




J13 Length 15.6cm, Tip Wide 1.0mm, Bayonet, non-stick epistaxis, Serrated Tip

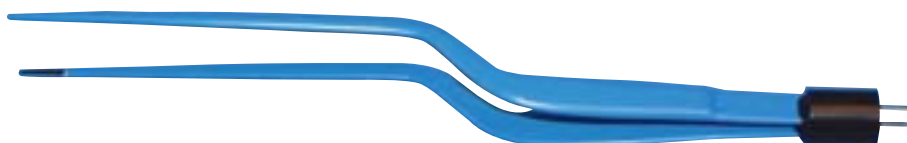

J17 Length 17.9cm, Tip Wide 1.0mm, Bayonet, external tip insulated

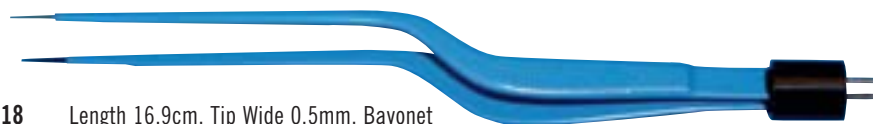

J18 Length 16.9cm, Tip Wide 0.5mm, Bayonet

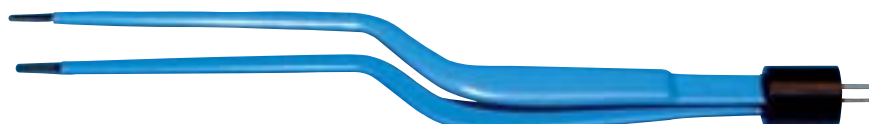

J20 Length 16.6cm, Tip Wide 2.0mm, Bayonet



J21 Length 19.7cm, Tip Wide 1.0mm, Bayonet




ellman® Bayonet Bipolar Forceps enhance access to the surgical environment for increased visibility during radiosurgical procedures.

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

RADIOURGICAL ACCESSORIES

ACe-Tip™ Non-Stick Forceps


Featuring Advanced Composition Alloy “Non-Stick”, Low Temperature

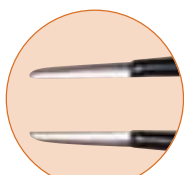

Magnified view

ACBF-010


Length (3.5"), 8.89cm, 0.5mm, Straight

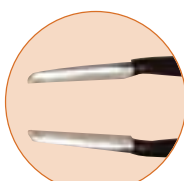

Magnified view

ACBF-011


Length (3.5"), 8.89cm, 0.5mm, Bent

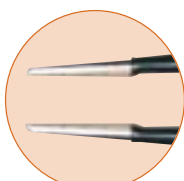

Magnified view

ACBF-012


Length (4.5"), 11.43cm, 0.5mm, Straight

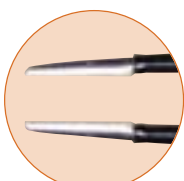

Magnified view

ACBF-013

Length (5.5"), 14.00cm, 0.5mm, Straight

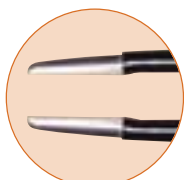

Magnified view

ACBF-014

Length (6.25"), 15.88cm, 1.0mm, Straight



ACe-Tip™ Technology
Advantages:*

- Minimal adherence of tissue
- Low absorption
- Tissue preservation




For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

*Data on file.
FDA Submission Packet
510 (k) K071343

RADIOSURGICAL ACCESSORIES

ACe-Tip™ Non-Stick Forceps


Featuring Advanced Composition Alloy “Non-Stick”, Low Temperature

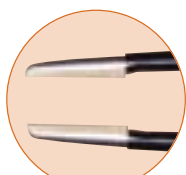

Magnified view




ACBF-015 Length (7.0”), 17.18cm, 0.7mm, Straight

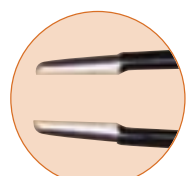

Magnified view




ACBF-016 Length (7.0”), 17.18cm, 1.2mm, Straight

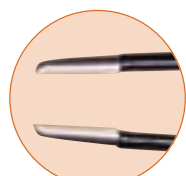

Magnified view




ACBF-017 Length (7.0”), 17.18cm, 1.5mm, Straight

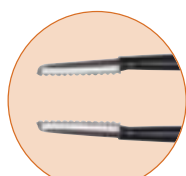

Magnified view



ACBF-021 Length (6.5”), 16.51cm, 1.0mm, Bayonet, Serrated Tip

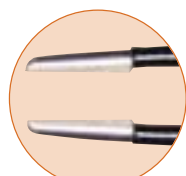

Magnified view



ACBF-022 Length (7.25”), 18.52cm, 1.0mm, Bayonet



ACe-Tip™ Non-Stick Bipolar
Forceps clinical photos

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

RADIOSURGICAL ACCESSORIES

RF SurgiPen™



SP010 # 4 Curette, 4mm
Sterile/Single-Use, (15 pcs)



SP011 # 7 Curette, 7mm
Sterile/Single-Use, (15 pcs)



SP012 Ball
Sterile/Single-Use, (15 pcs)



SP013 Loop
Sterile/Single-Use, (15 pcs)



SP014 Needle
Sterile/Single-Use, (15 pcs)



SP015 # 15 Scalpel
Sterile/Single-Use, (15 pcs)



SP016 Round Blade
Sterile/Single-Use, (15 pcs)



SP017 Vari-Tip™
Sterile/Single-Use, (15 pcs)

Quick-Connect™ Cable

IEC-SP020

RF SurgiPen™ to be used with the Quick-Connect™ Cable.
See page 11 for more information.

The Patented RF SurgiPen™ can be energized or simply used as a standard instrument.


Most surgical specialties will benefit from the use of this line of unique instruments due to their precision and versatility.



SurgiPen's™ ergonomic handle design is composed of lightweight material for fatigue-free manipulation.


For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

RF ELECTRODES

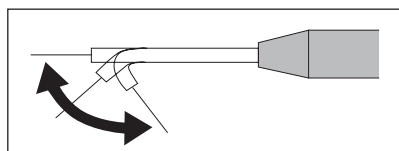
Electrodes



	(3/32") 2.4mm
	(1/16") 1.6mm

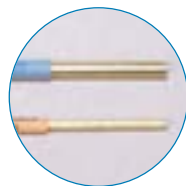

Proximal Shaft Diameter

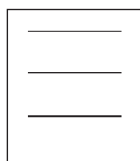

Super Fine - Used when minimal tissue alteration is desired (Diam. 0.005")

Fine - Used when less tissue alteration is desired (Diam. 0.007")

Regular - Used to produce more coagulation (Diam. 0.009")

ABOUT

ellman® Electrodes:

The insulated shaft of most electrodes can be bent or contoured for all surgical procedures.

(Note: Empire microIncision™ Sterile Needles are not bendable)

 **Yellow** electrode tips are Sterile/Single-Use disposables

 **Blue** electrode tips are re-usable

 **Orange** electrode tips are re-usable

 **Purple** electrode tips are re-usable

 **White** electrode tips are a mixture of re-usable and disposable.

See product catalog for details

INFORMATION

Regarding Ordering Electrodes:

Code To Electrode Suffixes:

B- Bendable

D- Disposable (Sterile/Single-Use)

Code To Electrode Prefixes:

T- Package of (2 pcs)

X- Package of (5 pcs)

Most electrodes are available in 1/16" or 3/32" shaft sizes.

Use "Q" as a prefix when ordering 3/32" shaft size, i.e. QTA8B

Longer electrode shaft sizes are available.

Contact **ellman's** sales representatives for more details.

Proximal shaft diameter is available in 3/32" or 1/16" (2.4mm or 1.6mm) to accommodate standard handpiece connection.

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

RF ELECTRODES

Sterile/Single-Use

Medical Electrode Set

DSEP40

Sterile/Single-Use

Assortment of 8 styles, Sterile/Single-Use










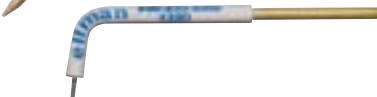
RF electrodes, available in 1/16" Shaft Only (40 pcs)

Includes 5 of each electrode:

A2D, A8D, B1D, B2D, C3D, D8D, E3D, F1D

Sterile/Single-Use Electrodes

Available in 1/16" Shaft Only (25 pcs per pack)

	A2D	Fine Wire Electrode 3/8" Length
	A3D/OS	Needle Electrode 3/4" Length
	A8D	Vari-Tip™ Wire Electrode Adjustable cutting wire
	B1D	Round Loop Electrode 1/4" Dia.
	B2D	Round Loop Electrode 3/8" Dia.
	C3D	Diamond Loop Electrode 3/16" Dia.
	D8D	Ball Electrode 2.3 mm
	E3D	Scalpel Electrode 1" Length
	E4D	Scalpel Electrode 1/2" Length
	F1D	Broad Needle Electrode
	133D	Pin Electrode 1/8" Length

The **ellman®** Sterile/Single-Use RF electrode eliminates the risk of cross contamination as well as ensuring consistent performance and patient confidence.

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com


RF ELECTRODES

Empire® microIncision™ Needles

Empire microIncision™ Needle Electrodes

Sterile/Single-Use, available in 1/16" or 3/32" (2 per pack)

Use "Q" as a prefix when ordering 3/32" shaft size

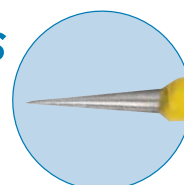

Magnified view
for all Empire
microIncision™



TEE301
Length 2cm, straight



TEE302
Length 3cm, straight



TEE303
Length 4cm, straight



TEE304
90° angle



TEE305
45° angle



TEE306
Length 3cm, 3mm 45° angle



TEE307
Length 2cm, 3mm 45° angle



TEE308 0.5cm needle,
Length (4"), 10.10cm, straight sleeve



TEE309 0.5cm needle,
Length (5"), 12.70cm, straight sleeve



TEE310 0.5cm needle,
Length (6"), 15.24cm, straight sleeve



TEE311 0.5cm needle,
Length (7"), 17.78cm, straight sleeve



TEE312 0.5cm needle,
Length (4"), 3mm 45° angle, straight sleeve



TEE313 0.5cm needle,
Length (5"), 3mm 45° angle, straight sleeve



TEE314 0.5cm needle,
Length (6"), 3mm 45° angle, straight sleeve



TEE315 0.5cm needle,
Length (7"), 3mm 45° angle, straight sleeve

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

The ellman® Empire®

microIncision™ Needle

Electrodes deliver focused

High Frequency Radiowave

energy.

Empire® microIncision™

Needle Electrodes are used for

precise cutting and coagulation

during radiosurgical procedures.

Proximal shaft diameter is

available in 3/32" or 1/16"

(2.4mm or 1.6mm) to

accommodate standard

handpiece connection.

Empire® microIncision™

Needle Electrodes are supplied

in individual sterile packages

of two, are Single-Use and

should not be re-sterilized.


RF ELECTRODES

Empire[®] microIncision[™] Needles

Empire microIncision[™] Needle Electrodes

Sterile/Single-Use, available in 1/16" or 3/32" (2 per pack)


Use "Q" as a prefix when ordering 3/32" shaft size

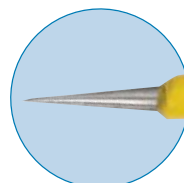

Magnified view
for all Empire
microIncision[™]



TEE316 0.5cm needle,
Length (4"), 10.16cm, 45° angle, straight sleeve

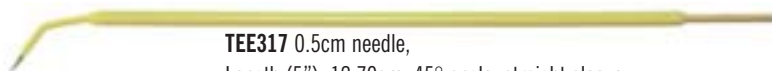

TEE317 0.5cm needle,
Length (5"), 12.70cm, 45° angle, straight sleeve



TEE318 0.5cm needle,
Length (6"), 15.24cm, 45° angle, straight sleeve




TEE319 0.5cm needle,
Length (7"), 17.78cm, 45° angle, straight sleeve


TEE320 0.5cm needle,
Length (4"), 10.16cm, curved sleeve



TEE321 0.5cm needle,
Length (5"), 12.70cm, curved sleeve

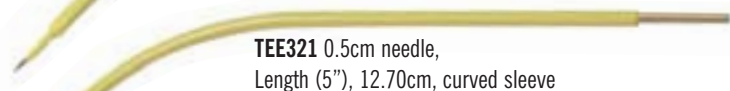

TEE322 0.5cm needle,
Length (6"), 15.24cm, curved sleeve

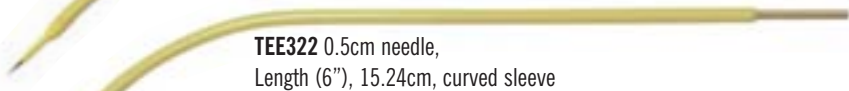

TEE323 0.5cm needle,
Length (7"), 17.78cm, curved sleeve



TEE324 Bayonet



TEE325 0.5cm needle,
Length (2"), 5.08cm, curved sleeve



TEE326 0.5cm needle,
Length (1"), 2.54cm,
40 degree curved sleeve



TEE327 0.5cm needle,
Length (3/4"), 1.91cm, straight



TEE328 0.5cm needle,
Length (1"), 2.54cm,
30 degree curved sleeve



Upper Lid Blepharoplasty



For orders and/or pricing information, please call ellman[®] at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

RF REUSABLE ELECTRODES

Reusable Electrodes

Needle Electrodes*

Reusable

1/16" shaft size or 3/32" (2 pcs)

Use "Q" as a prefix when ordering 3/32" shaft size



- | | |
|--------------|--|
| TA1B | Length 3/8", Super Fine Wire |
| TA2B | Length 3/8", Fine Wire |
| TA9B | Length 3/8", Regular Wire |
| TA3B | Length 3/8", Regular Needle |
| TF1B | Length 1/2", Broad Needle |
| TA8B | Vari-Tip™ Microdissection™ Needle (comes with additional wire and has an adjustable depth) |
| T133B | Pin Electrode 1/8" Length |

Round Loop Electrodes*

Reusable

1/16" shaft size or 3/32" (2 pcs)

Use "Q" as a prefix when ordering 3/32" shaft size



- | | |
|--------------|--------------------------------------|
| TBBB | 5/32" Dia. Round Loop, Fine Wire |
| TB1B | 1/4" Dia. Round Loop, Fine Wire |
| TB4B | 1/4" Dia. Round Loop, Regular Wire |
| TB9B | 19/64" Dia. Round Loop, Fine Wire |
| TB10B | 19/64" Dia. Round Loop, Regular Wire |
| TB2B | 3/8" Dia. Round Loop, Fine Wire |
| TB5B | 3/8" Dia. Round Loop, Regular Wire |
| TB3B | 1/2" Dia. Round Loop, Fine Wire |
| TB6B | 1/2" Dia. Round Loop, Regular Wire |

**Electrodes are displayed in actual size.*

Needle Electrodes are commonly used for incisions and excisions in precise surgical procedures.

Choice of Wire Gauge

- Super Fine - Used when minimal tissue alteration is desired.
- Fine - Used when less tissue alteration is desired.
- Regular - Used to produce more coagulation.

Round Loop Electrodes are commonly used for shaving lesions, contouring tissue and obtaining biopsy specimens.

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

RF REUSABLE ELECTRODES

Reusable Electrodes

Diamond Loop Electrodes*

Reusable

1/16" shaft size or 3/32" (2 pcs)

Use "Q" as a prefix when ordering 3/32" shaft size



TC1B 1/8" Diamond Loop, Fine Wire

TC4B 1/8" Diamond Loop, Regular Wire



TC3B 3/16" Diamond Loop, Fine Wire

TC6B 3/16" Diamond Loop, Regular Wire



TC17B 1/4" Diamond Loop, Fine Wire

TC18B 1/4" Diamond Loop, Regular Wire



TC19B 5/16" Diamond Loop, Fine Wire

TC20B 5/16" Diamond Loop, Regular Wire



TC21B 7/16" Diamond Loop, Fine Wire

TC22B 7/16" Diamond Loop, Regular Wire

Triangle Loop Electrodes*

Reusable

1/16" shaft size or 3/32" (2 pcs)

Use "Q" as a prefix when ordering 3/32" shaft size



TC9B 5/32" Triangle Loop, Fine Wire

TC10B 5/32" Triangle Loop, Regular Wire



TC7B 1/4" Triangle Loop, Fine Wire

TC8B 1/4" Triangle Loop, Regular Wire



TC11B 19/64" Triangle Loop, Fine Wire

TC12B 19/64" Triangle Loop, Regular Wire



TC13B 3/8" Triangle Loop, Fine Wire

TC14B 3/8" Triangle Loop, Regular Wire



TC15B 1/2" Triangle Loop, Fine Wire

TC16B 1/2" Triangle Loop, Regular Wire

**Electrodes are displayed in actual size.*

Diamond Loop Electrodes

are commonly used for

one-stroke elliptical incisions.

Triangle Loop Electrodes

are commonly used for

shaving lesions and

contouring tissue.

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

RF REUSABLE ELECTRODES

Reusable Electrodes

Ball Electrodes*

Reusable

1/16" shaft size or 3/32" (2 pcs)

Use "Q" as a prefix when ordering 3/32" shaft size



TD3B	2mm Ball
TD8B	2.36mm Ball
TD9B	3.17mm Ball
TD10B	4.34mm Ball
TD11B	4.75mm Ball

Oval Loop Electrodes*

Reusable

1/16" shaft size or 3/32" (2 pcs)

Use "Q" as a prefix when ordering 3/32" shaft size



TP1B	5/32" Oval Loop, Fine Wire
TP2B	5/32" Oval Loop, Regular Wire
TP3B	1/4" Oval Loop, Fine Wire
TP4B	1/4" Oval Loop, Regular Wire
TP5B	19/64" Oval Loop, Fine Wire
TP6B	19/64" Oval Loop, Regular Wire
TP7B	3/8" Oval Loop, Fine Wire
TP8B	3/8" Oval Loop, Regular Wire
TP9B	1/2" Oval Loop, Fine Wire
TP10B	1/2" Oval Loop, Regular Wire

**Electrodes are displayed in actual size.*

Ball Electrodes are commonly used for coagulation.

Oval Loop Electrodes are commonly used for shaving lesions and contouring tissue.

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

RF REUSABLE ELECTRODES

Reusable Electrodes

Scalpel Electrode*

Reusable

1/16" shaft size or 3/32" (2 pcs)

Use "Q" as a prefix when ordering 3/32" shaft size



TE1B 3/16" Scalpel, 45° Angle

TE2B 3/8" Scalpel

TE3B 1" Scalpel

Empire® MicroTip Electrodes*

Reusable


2mm exposed tip (2 pcs, one of each)



TEE227 Empire MicroTip Needles
one .010 wire, one .020 wire
Insulated Coated Needle Set

Insulated Needle Electrodes™*

Reusable, 1/2mm exposed tip

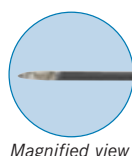

Magnified view



D7 Insulated Coated Needle Set
(3 pcs)

D6A Fine Insulated Coated Needle
.004 Super Fine, White, (1 pc)

D6B Fine Insulated Coated Needle
.007 Fine, Purple, (1 pc)

D6C Fine Insulated Coated Needle
.009 Regular, Orange, (1 pc)

Rhinophyma Electrode Set

EE234 Set of 3

Reusable



EE233 19/64" Triangle Loop, (1 pc)

EE232 1/8" Triangle Loop, (1 pc)


EE231 1/2" Triangle Loop, (1 pc)

**Electrodes are displayed in actual size.*

Scalpel Electrodes are commonly used for incisions and excision of tissue. The blade shape provides broader coagulation than a needle electrode.

Empire® MicroTip Needles are tapered and have a pinpoint design which produces micro-smooth cutting and coagulation while creating less tissue alteration.

Insulated Needle Electrodes are used for the treatment of Telangiectasia and for epilation.

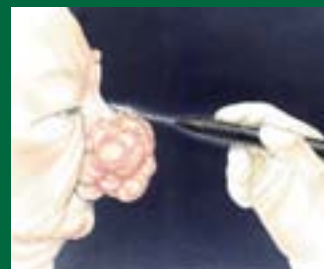

Rhinophyma


For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

RF REUSABLE ELECTRODES

Reusable Electrodes

Javate Endoscopic DCR Electrodes*


R6 Javate Endoscopic DCR Set of 7 Assorted Electrodes

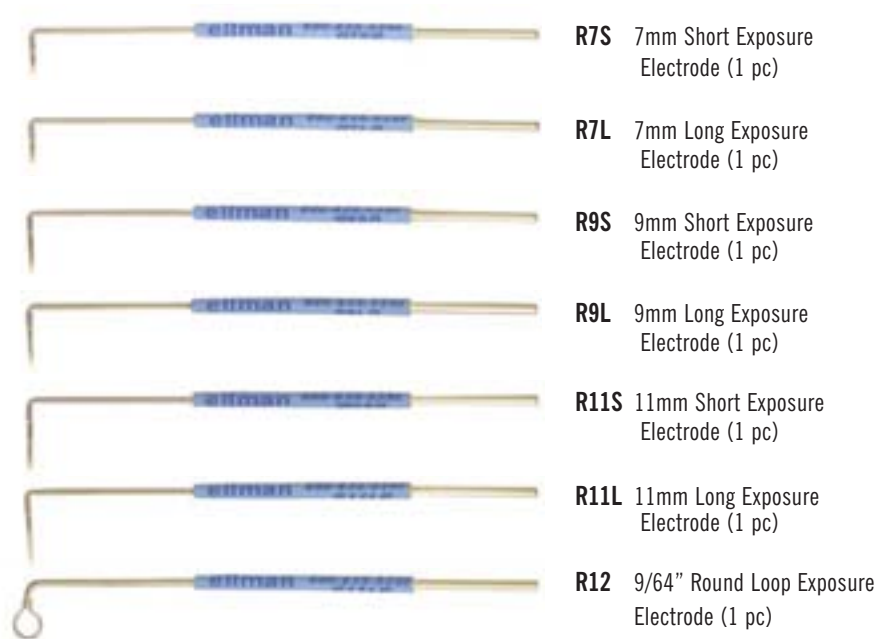

**Electrodes are displayed in actual size.*

Radiosurgical Mucotome Electrodes*

H140 Radiosurgical Mucotome Electrode Set of 1/16" Shaft Reusable
Set of 3 electrodes, includes H141, H142 & H143
Used for excision of mucosal strip

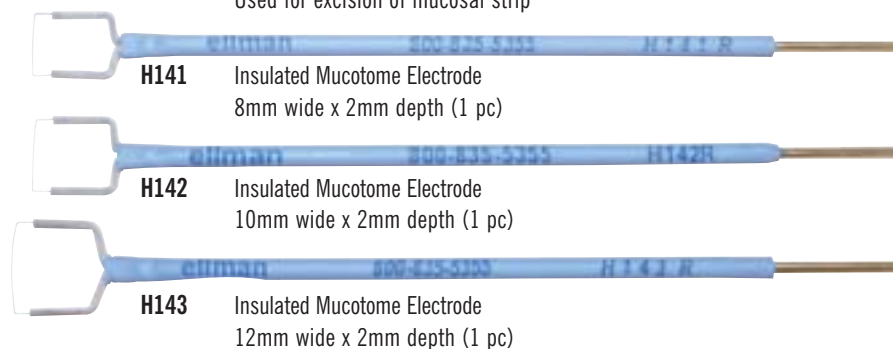

**Electrodes are displayed in actual size.*

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com


RF REUSABLE ELECTRODES*

ACe-Tip™ Advanced Composition Alloy

TNAEE286 HLTZ Tonsillotomy Set of 3 ACe-Tip™ electrodes
Includes one each of: TNAEE284, TNAEE285, TNAEE287



TNAEE284 HLTZ Tonsillotomy, ACe-Tip™,
15mm Loop, .010 wire, 1/16" shaft (2 pcs)

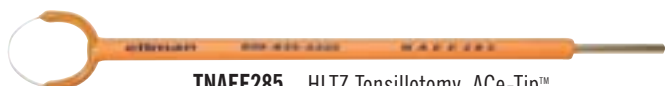

TNAEE285 HLTZ Tonsillotomy, ACe-Tip™,
12mm Loop, .010 wire, 1/16" shaft (2 pcs)



TNAEE287 1cm Needle, 3" Curved Sleeve ACe-Tip™, 1/16" shaft (2 pcs)



TNAA9B 3/8" Length
ACe-Tip™ Regular Wire 1/16" shaft (2 pcs)



TNAD3B 2mm Coagulating Ball
ACe-Tip™ 1/16" shaft (2 pcs)



TNATB3B 3mm Tapered Ball
ACe-Tip™ 1/16" shaft (2 pcs)



TNAD13B 5mm Tapered Ball
ACe-Tip™ 1/16" shaft (2 pcs)



TNAB1B 1/4" Round Loop
ACe-Tip™ Fine Wire 1/16" shaft (2 pcs)



TNAB4B 1/4" Round Loop
ACe-Tip™ Regular Wire 1/16" shaft (2 pcs)



TNAB2B 3/8" Round Loop
ACe-Tip™ Fine Wire 1/16" shaft (2 pcs)



TNAB5B 3/8" Round Loop
ACe-Tip™ Regular Wire 1/16" shaft (2 pcs)



TNAA8B Vari-Tip™
ACe-Tip™ 1/16" shaft (2 pcs)



TNAEE305 Empire microIncision™
Needle ACe-Tip™ 1/16" shaft (2 pcs)

ACe-Tip™ Technology

Advantages:*

- Minimal adherence of tissue
- Low absorption
- Tissue preservation

The **ellman® ACe-Tip™** electrodes

are designed to reduce resistance/impedance which causes less undesired lateral heat compared to tungsten steel electrodes.

A 45% reduction in the zone of hyalinization has been demonstrated.†

ACe-Tip™ electrodes glide

smoothly through tissues and are excellent for sculpting and tissue shaping.

*Data on file.

FDA Submission Packet
510 (k) K071343

†Porcine skin.

For orders and/or pricing information, please call **ellman®** at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

RF REUSABLE ELECTRODES

Reusable Electrode Sets

Rodin® Microsurgery ACe-Tip™ Electrode Set

RACE7

Reusable

1.5" shaft length, 1/16" dia., Bendable Assortment (7 pcs)

Includes one each of:

REE305, RA2B, RA8B, RD3B, RBBB, RC3B, RB1B



Rodin® Microsurgery Handpiece

Standard Electrode Set

S10B

Reusable

1/16" shaft size or 3/32" Use "Q" as a prefix when ordering 3/32" shaft size

Assortment (7 pcs) Includes one each of: C7B, B1B, C3B, A3B, D8B, F1B, A8B



Assorted Needle Electrode Set

XA

Includes one each of: A1, A2, A9, A3, F1 (5 pcs)



Assorted Round Loop Electrode Set

XB


Fine Wire

Includes one each of: B3, B2, B9, B1, BB

XBRB

Regular Wire

Includes one each of: B6, B5, B10, B4, BB (5 pcs)

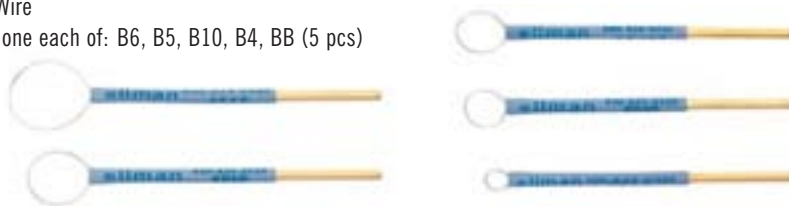

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

Reusable Electrodes Sets

Assorted Diamond Loop Electrode Set

XCFB

Fine Wire

Includes one each of: C21, C19, C17, C3, C1

XCRB Regular Wire

Includes one each of: C22, C20, C18, C6, C4
(5 pcs)



Assorted Ball Electrode Set

XD

Includes one each of: D11, D10, D9, D8, D3
(5 pcs)



Assorted Oval Loop Electrode Set

XPFB

Fine Wire

Includes one each of: P9, P7, P5, P3, P1

XPRB

Regular Wire

Includes one each of: P10, P8, P6, P4, P2
(5 pcs)



Assorted Triangle Loops Electrode Set

XTFB

Fine Wire

Includes one each of: C15, C13, C11, C9, C7

XTRB

Regular Wire

Includes one each of: C16, C14, C12, C10, C8
(5 pcs)



For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

RF REUSABLE ELECTRODES

Electrodes Accessories

Electrode Protector Stand™

H66

Electrode Protector Stand™

For regular 2" standard electrodes

Holds 12 electrodes, 6 1/16" x 6 3/32"

(1 pc)

H66A

Electrode Protector Stand™

For large 4" long electrodes

Holds 12 electrodes, 6 1/16" x 6 3/32"

(1 pc)



The **ellman®** Electrode Protector Stand™ protects your electrodes from damage and breakage during the sterilization process.

The **ellman®** Electrode Protector Stand™ will stand on any flat surface for easy identification and access to each electrode.

Electrode Block™

H1

Electrode Block™

Holds 18, 1/16" Electrodes

(1 pc)

H1A

Electrode Block™

Holds 9, 1/16" Electrodes and

9, 3/32" Electrodes

(1 pc)

QH1

Electrode Block™

Holds 18, 3/32" Electrodes

(1 pc)



For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

EDUCATIONAL MATERIALS

Educational Accessories

Compact Discs
Oculoplastic Techniques
by Jay Justin Older, MD
VCD/OP31
2 Disk Set



At **ellman®** we are committed to the ongoing education of clinicians as well as their patients.

Patient Education Brochures



Cosmetic Eye Surgery
LIT/OP1A
(50 pcs)



Facial Veins
LIT/DE1
(50 pcs)



Skin Tightening*
LIT/DE2
(50 pcs)




Facial Lesions
LIT/CO1
(50 pcs)



Treatment of Snoring
LIT/ENT1
(50 pcs)

Textbooks

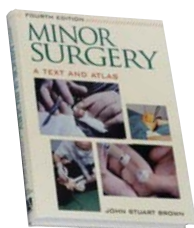

Minor Surgery
H202A
A Text & Atlas
4th Edition,
by John Stuart Brown



Procedures for Primary Care
2nd Edition
H208
by Drs. John L. Pfenninger and Grant C. Fowler



Eyelid Tumors
H207
by Jay Justin Older, MD



Oculoplastic Surgery with Radiofrequency
H205/BL
by Drs. G. Aimino, G. Davi, and M. Santella

***FDA 510(k) pending. Skin Tightening product is only available to physicians outside the United States.**

For orders and/or pricing information, please call **ellman®** at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

SMOKE EVACUATION SYSTEM

Device/Accessories



Patented Surg-e-Vac™

SV110 Surg-e-Vac™ (for 110v countries)

SV220EU Surg-e-Vac™ (for 220v countries)

SurgiCart™

H156

SurgiCart™ III with extended base, instrument tray with gooseneck arm, adjustable height, hospital locking wheels, rolls easily from room to room



Vac Assist Stand

VV133

The Surg-e-Vac™ Smoke Evacuation System has been designed with a high suction, high flow rate centrifugal action pump. The ultra-quiet motor is used to draw the surgical smoke from the surgical site through the vacuum tubing and into the Quad Filter™ where the surgical smoke is processed by a series of filters.

ellman® brings you a strong smoke suction system, the Surg-e-Vac™.

The unique compact multi-filtration technology has been designed to trap macro and micro airborne particles.

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

SMOKE EVACUATION SYSTEM

Accessories

Surg-e-Vac™ Accessories

The accessories below are used only for Surg-e-Vac™

Surg-e-Vac™
Quadration
Filter
SVF1 (3 pcs)
SVXF1 (12 pcs)



Pneumatic
Footswitch
SVFSC
(10-1/2 Ft. Cable)
(1 pc)



For both Surg-e-Vac™ and Vapor-Vac II™ Accessories

Plastic Hose
VV148 or VV127
Reusable
7/8" diameter
(1 pc)



Speculum Tubing
Reducing Package
VV130 includes
(2 pcs) VV128
(2 pcs) VV129



Saf-T-Wand
VV140
(10 pcs)



Speculum Tubing
Reducer Fitting
VV129
(2 pcs)



Speculum Tubing
VV128
1/4" Autoclavable
(2 pcs)




External Tube
VV121
(3 pcs)



Vapor-Vac II™ Accessories

These accessories below are used only for Vapor-Vac II™

Charcoal
Filter
VC123
(2 pcs)

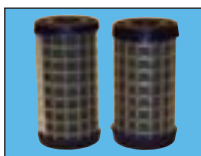

External
Filter
VF124
(10 pcs)




Saf-T-Wand
Filter Set
VV146
(25 pcs)



Disposable
External Filter
VV139
(25 pcs)



Vapor-Vac™
Wand Clip
VV131
(1 pc)

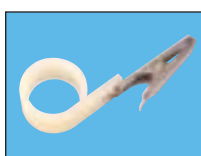

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

MULTI-SPECIALTY ACCESSORIES

Accessories

Medi-Frig™

Portable Cryosurgery Skin Refrigerant

H11B

1 280g. can, 10 Straw Applicators

XH11B

6 cans, 60 Straw Applicators

H11BCS


25 cans, 250 Straw Applicators



Wipe-Its™

H19

Pressure sensitive sandpaper pads
used to clean radiosurgical electrode tips
(4 pcs)

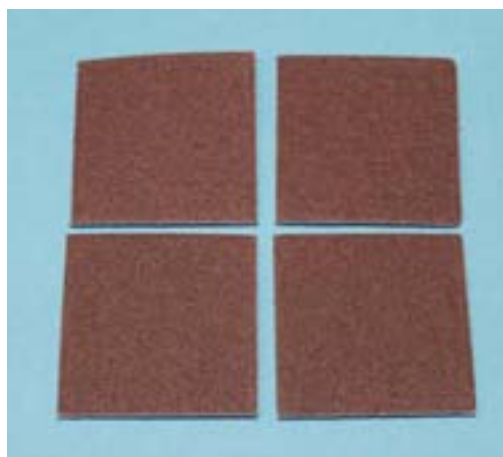

Medi-Frig™ is a portable, hand-held, cryosurgical spray unit containing tetrafluoroethane (-54°C).

The **ellman**® cryosurgical spray is safe, storable, easily transportable and cost-effective compared to traditional liquid nitrogen. An added advantage is that it is nonflammable and can be used in conjunction with radiosurgical, electrosurgical and laser procedures.

Wipe-Its™ are easy to use, for effortless removal of tissue from radiosurgery tips.

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

MULTI-SPECIALTY ACCESSORIES

Otorhinolaryngology



3" Length Handpiece with Patented Locking Chuck

IEC-XLP1

Reusable

For Surgitron® 4.0 Dual RF™ and Radiolase® II Devices

3" Length, Accepts 1/16" and 3/32" shaft, procedures for Laparoscopy and Micro-Larynx (2 pcs)

XLP1

Reusable

For the EMC Device

3" Length, Accepts 1/16" and 3/32" shaft, for

Laparoscopy and Micro-Larynx procedures (2 pcs)




Micro-Larynx™ RF Probes

IEC-EE275A

Reusable

For Surgitron® 4.0 Dual RF™ and Radiolase® II Devices

1/16" shaft, 240mm long, (3 pcs)

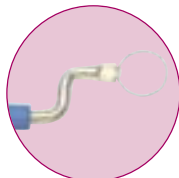


Magnified view

Loop

EE276

1/16" shaft, 240mm long
(1 pc)

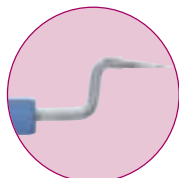


Magnified view

Needle

EE277

1/16" shaft, 240mm long
(1 pc)

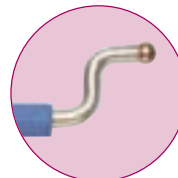

Magnified view

Ball

EE278A

1/16" shaft, 240mm long
(1 pc)

Patented Micro Fibre® Electrode

XMF-A

Sterile/Single-Use

Monopolar Wire

600mm Length (1mm O.D.)

for endoscopic procedures,

Attaches to Quick-Connect™ Cable (3 pcs)



Magnified view



The Patented Micro-Larynx™ Probe allows quick and efficient coagulation, and can be used for Granulomas, Papillomatosis, Cyst Biopsies, Vocal Cord Polyps, Posterior Cordectomy, Glottic and Sub-glottic Angiomas, Laryngeal Tumor and Laryngeal Biopsies.

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

MULTI-SPECIALTY ACCESSORIES

Otorhinolaryngology

OTO-RAM™ Radiofrequency Assisted Myringotomy Set

IEC-EE600

Reusable


Two 4" tubes, 1.5mm and 2.0mm and one Myringotomy Handpiece

IEC-H604

Reusable

RF Assisted Myringotomy Handpiece (1 pc)

The OTO-RAM™ allows for an effective treatment where office or clinical procedures take place in only a few minutes. Costs for OR Surgery and extensive use of antibiotics and general anesthesia are greatly reduced. The OTO-RAM™ is cost effective compared to laser assisted myringotomy.

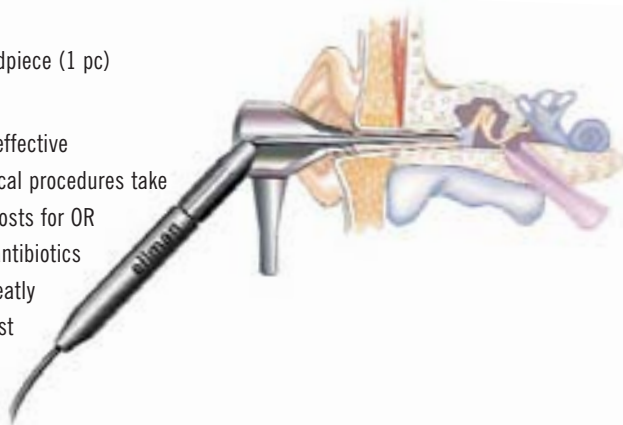

Patented OTO-RAM™ Probes

EE601

Reusable

1.5mm Orifice


(1 pc)

EE602

Reusable

2.0mm Orifice

(1 pc)

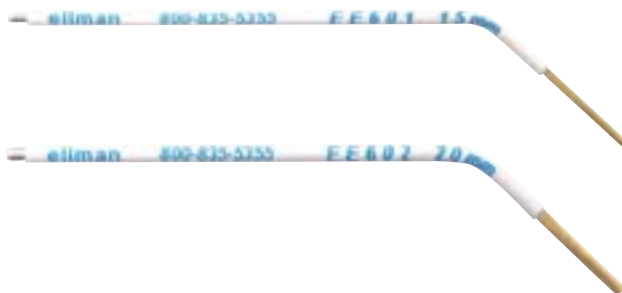

Patented RF Myringotomy Electrode™

TEE230

Reusable

Insulated 1cm and exposed 1cm

(2 pcs)




Quick-Connect™ Cable

Reusable

Proximal cable for Micro Fibre® System.

See page 11 for ordering information

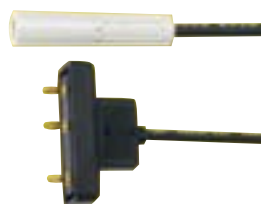

2mm hole



2mm hole
Tympanic membrane



OTO-Ram
Probe Tip

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

MULTI-SPECIALTY ACCESSORIES

Otorhinolaryngology

Clear-Vu® Bayonet Electrodes

Sterile/Single-Use, (13cm Long Shaft), (1 pc)

1/16" shaft size, 3/32" shaft size available, order with a "Q" prefix



Bayonet RF Turbinate Volume Reduction Electrode™

Reusable

Bipolar Turbinate Probe and cord (1 pc)

Tips Insulated 1cm and Exposed 1cm

IEC-I6D

For Surgitron® 4.0 Dual RF™ and Radiolase® II Devices

I6/JXE

For EMC Device



Magnified view

The bayonet design of Clear-Vu® Electrodes provides access to the surgical environment for precise tissue interaction while maintaining a 'clear view' during the procedure.

The assorted tip configurations allows for precise surgical treatment.



The Bayonet Turbinate Electrode and Bipolar Handpiece design includes a pencil grip for easy handling with a bayonet needle for easy intranasal visualization during surgery.

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

MULTI-SPECIALTY ACCESSORIES

Otorhinolaryngology

Patented RAUP Electrode™

H81A1


Sterile/Single-Use
RAUP Radiosurgically Assisted
Uvulopalatoplasty
(2 pcs)



Sub-Mucosal VTR (UPP) Electrode

TEE229A

Sterile/Single-Use
For volumetric reduction
of redundant tissue
(2 pcs)

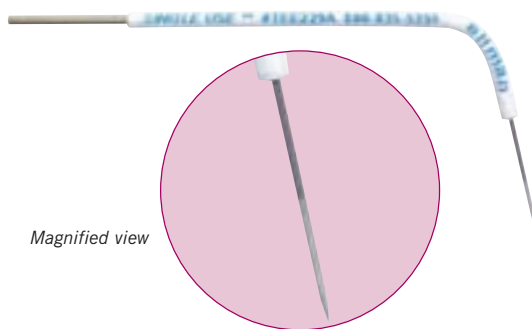

Magnified view

Patented Tongue Base Reduction Electrode™

TEE605

Sterile/Single-Use,
For volumetric reduction
of redundant tissue
(2 pcs)



Magnified view

Adenoid RF Electrode™

TEE228


Reusable
1/4" Triangle, 4" Length, Teflon Insulated
(2 pcs)



HLTZ Tonsillotomy Electrode™

TEE284

Reusable
.010 Wire, 15mm Loop
(2 pcs)

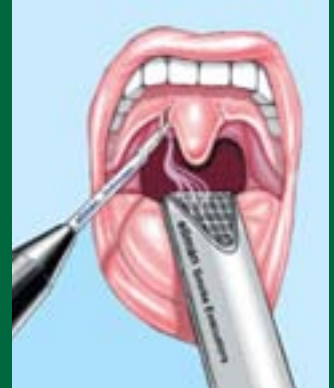

Uvulopalatoplasty (RAUP)



Submucosal Stiffening



Tongue Base Procedure

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

MULTI-SPECIALTY ACCESSORIES

Otorhinolaryngology



Patented RF Suction Adenoid Probes

IEC-H85

Reusable

Set of 3 Probes, Quick-Connect™ Cable, and Cleaning Brush



IEC-H86 Small Probe
9mm with Quick-Connect™ Cable, (1 pc)



IEC-H87 Medium Probe
15mm with Quick-Connect™ Cable, (1 pc)



IEC-H88 Large Probe
19mm with Quick-Connect™ Cable, (1 pc)


Cleaning Brush for Suction Probe

H89

(3 pcs)



Magnified view

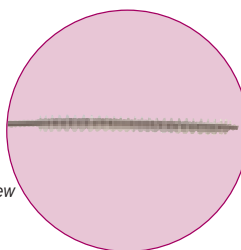

Patented Suction Tonsillectomy Probe™ and Quick-Connect™ Cable

IEC-H83/SP

Reusable

For Surgitron® 4.0 Dual RF™ and Radiolase® II Devices

H83/SP

For EMC Device



Patented RF Bipolar Tonsil Reduction Probe™ and Quick-Connect™ Cable

IEC-H84

Reusable

For Surgitron® 4.0 Dual RF™ and Radiolase® II Devices

H84

For EMC Device

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

MULTI-SPECIALTY ACCESSORIES

Otorhinolaryngology

Patented RF Suction Coagulator™

HD41

Sterile/Single-Use

For use with Quick-Connect™ Cable, (12 pcs)
(see page 11 for Quick-Connect™ Cable)



Magnified view

Patented RF Angled Suction Coagulator™

HD42

Sterile/Single-Use

For use with Quick-Connect™ Cable, (12 pcs)
(see page 11 for Quick-Connect™ Cable)



Magnified view



Patented Suction Coagulation Kit

Reusable

IEC-H40

Shipped with three 3" and three 6" tubes
For Surgitron® 4.0 Dual RF™ and Radiolase® II
Devices. Includes 1/16", 3/32" and 1/8" Orifice sizes

H40

Shipped with three 3" and three 6" tubes
For EMC Device. Includes 1/16", 3/32"
and 1/8" Orifice sizes



Stylus Round Wire

H44 Suction Tube 6" Long, 1/16" Orifice

H45 Suction Tube 6" Long, 3/32" Orifice

H46 Suction Tube 6" Long, 1/8" Orifice

H47 Suction Tube 3-3/4" Long, 1/16" Orifice

H48 Suction Tube 3-3/4" Long, 3/32" Orifice

H49 Suction Tube 3-3/4" Long, 1/8" Orifice

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

MULTI-SPECIALTY ACCESSORIES

Otorhinolaryngology

RF Assisted Tonsillotomy (RAT) Ball Electrodes™

EE264 Set of 3

Reusable


TEE261 2.0mm Ball 4" long (2 pcs)



Magnified view




TEE262 2.36mm Ball 4" long (2 pcs)

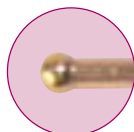

Magnified view



TEE263 3.17mm Ball 4" long (2 pcs)

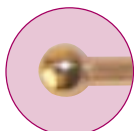

Magnified view




Cryptolysis Electrode

TE1022B

Reusable

(2 pcs)

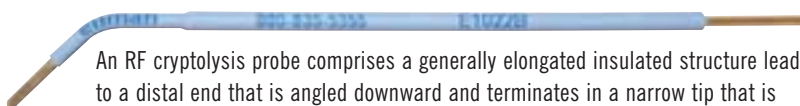

An RF cryptolysis probe comprises a generally elongated insulated structure leading to a distal end that is angled downward and terminates in a narrow tip that is exposed. The exposed tip is inserted into tonsillar crypts and activated.

Polyp Snare™ System for Nasal Polyps

IEC-H301

Reusable

For Surgitron® 4.0 Dual RF™ and Radiolase® II Devices
with Snare Cable and 3 wires


H301

For EMC Device
with Snare Cable and 3 wires

H7A Snare Wires
6 wires

H7B Snare Cable
For EMC Device (1 pc)

IEC-H7B Snare Cable
For Surgitron® 4.0 Dual RF™ and Radiolase® II
Devices (1 pc)

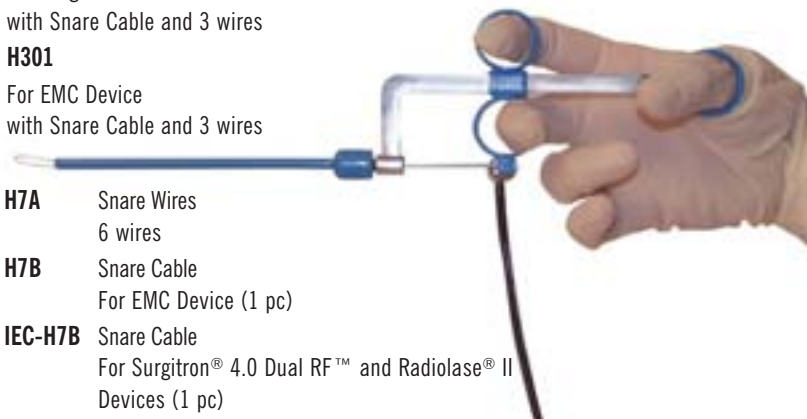

The Ball Tonsil Probe is an advancement over the traditional tonsillotomy. The RAT procedure permits controlled reduction of the tonsil without excessive thermal damage or blood loss.



Cryptolysis

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

MULTI-SPECIALTY ACCESSORIES

General Surgery/Proctology

Patented RF Laparoscopic Handpiece


IEC-XLP1

Reusable

3" Length Handpiece.

Accepts 1/16" and 3/32" shafts


For Surgitron® 4.0 Dual RF™ and Radiolase® II Devices
(1 pc)

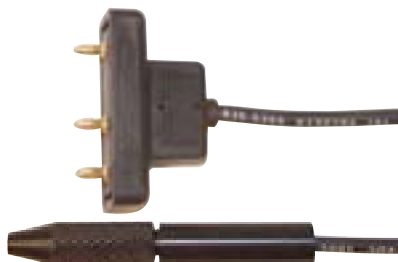

Patented RF Laparoscopic Electrodes™

Reusable

Lap-Trode Individual Electrodes, 37.5mm Length Electrodes, 1/16" Shaft

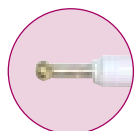

Magnified view

L3 2mm Ball Lap-Trode



Magnified view

L4 9.5mm Needle Lap-Trode




Magnified view

L5 2.3mm Hook Lap-Trode



Magnified view

L6 2.3mm Blade Lap-Trode

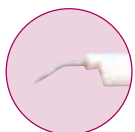

Magnified view

L8 Hook Needle Lap-Trode



Patented RF Laparoscopic Electrodes™ Set

L7

Reusable

Lap-Trode Set of 5 Assortment, (L3, L4, L5, L6, L7, L8) 37.5mm Length Electrodes, 1/16" Shaft

Patented Gupta RF Gun (Pravin Gupta Procto-Gun)

IEC-FSGUN

Reusable

Accepts 1/16" and 3/32" Shaft Electrodes (1 pc)

For Surgitron® 4.0 Dual RF™



RF Fingerswitch
Gun Handle

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

MULTI-SPECIALTY ACCESSORIES

Gynecology

LLETZ Electrodes Assortment Set

- W4** One of each in 5 styles, reusable (5 pcs)
W4D One of each in 5 styles, Sterile/Single-Use (5 pcs)

LLETZ Individual Electrodes

Available in 1/16" or 3/32"

Use "Q" as a prefix when ordering 3/32" shaft size

Use "D" as a suffix when ordering Sterile/Single-Use Versions of the LLETZ loop electrodes.

- W1** 15mm x 10mm Loop, reusable (1 pc)
W1D 15mm x 10mm Loop, Sterile/Single-Use (5 pcs)
- W2** 20mm x 10mm Loop, reusable (1 pc)
W2D 20mm x 10mm Loop, Sterile/Single-Use (5 pcs)
- W3** 25mm x 10mm Loop, reusable (1 pc)
W3D 25mm x 10mm Loop, Sterile/Single-Use (5 pcs)
- F2** Endo-Cervical Box 8mm x 11mm Box, reusable (1 pc)
F2D Endo-Cervical Box 8mm x 11mm Box, Sterile/Single-Use (5 pcs)
- W18** 5mm Ball, reusable (1 pc)
W18D 5mm Ball, Sterile/Single-Use (5 pcs)



Conization Electrode™

Reusable

- G1** Conization Electrode
G1/4 4" long, Rotated at 360° to provide a cone shaped specimen



Standard Electrode Set

S10B/4 4" shaft length, (7pcs) see page 33 for more information

External Genital Lesion RF Electrodes™

Reusable

- W15** Complete Set of Genital Lesion Electrode
 set of 4, includes W16, W17, W18 and W19 Electrodes

- W16** 90mm Length, Angled Box



- W17** 4mm Length, Box



- W18** 90mm Length, 5mm Ball



- W19** 90mm Length, 3/8" Needle



The **ellman®** LLETZ loop electrodes are used for the diagnosis and treatment of cervical intraepithelial neoplasia (CIN) and dysplasia.

For orders and/or pricing information, please call ellman® at: (800) 835-5355
 (516) 594-3333 • fax (516) 569-0054 • www.ellman.com

MULTI-SPECIALTY ACCESSORIES

Gynecology

LLETZ Bayonet Electrodes

- W35** LLETZ Bayonet Electrode Set (5 pcs)
Reusable. Use "D" as a suffix when ordering Sterile/Single-Use
- W35D** Sterile/Single-Use (5 pcs)
- W32** Reusable
Bayonet, 25mm Loop (1 pc)
- W32D** Sterile/Single-Use (5 pcs)
- W31** Reusable
Bayonet, 20mm Loop (1 pc)
- W31D** Sterile/Single-Use (5 pcs)
- W30** Reusable
Bayonet, 15mm Loop (1 pc)
- W30D** Sterile/Single-Use (5 pcs)
- W34** Reusable
Bayonette, 5mm Ball (1 pc)
- W34D** Sterile/Single-Use (5 pcs)
- W33** Reusable
Endo-Cervical Bayonet,
8mm Box (1 pc)
- W33D** Sterile/Single-Use (5 pcs)




Vaginal Speculums

- H55** Vaginal Speculum - Large, 5" Length
- H57** Vaginal Speculum - Standard, 4 1/2" Length
- H58** Vaginal Speculum (2 pcs - Large and Standard)



Vaginal Lateral Wall Retractor

H69

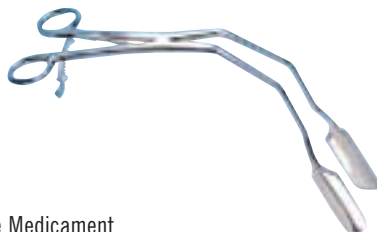

Monsel's Solution

H90

500 ml (Ferric Subsulfate) Post-Operative Medicament



Needle Extender/Measured Dose Syringe

H71




For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

MULTI-SPECIALTY ACCESSORIES

Plastic Surgery

Patented Javate/Khan RF Endoscopic Forehead Suction Probe™

Reusable

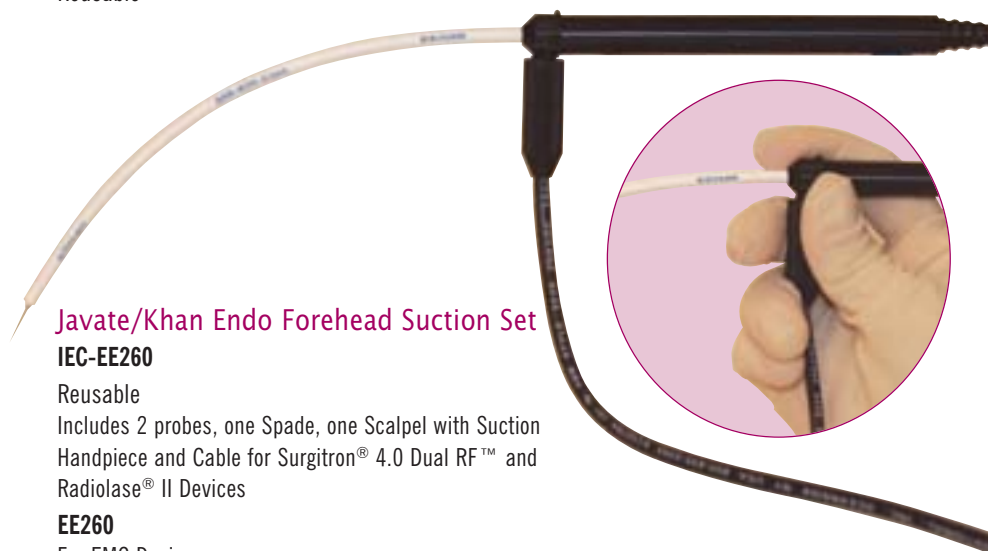

Javate/Khan Endo Forehead Suction Set

IEC-EE260

Reusable

Includes 2 probes, one Spade, one Scalpel with Suction Handpiece and Cable for Surgitron® 4.0 Dual RF™ and Radiolase® II Devices

EE260

For EMC Device




The Javate/Khan Endoscopic Forehead Suction Probe™ utilizes a combination of suction and radiofrequency cutting.

Javate/Khan Coated Flat Scalpel Suction Probe

EE257

Reusable

7/16" Tip Diameter (1 pc)

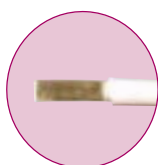

Magnified view



Javate/Khan Coated Spade Suction Probe

EE258

Reusable

7/16" Tip Diameter (1 pc)



Magnified view



For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

MULTI-SPECIALTY ACCESSORIES

Plastic Surgery

Patented Dr. Constantin Stan Monopolar Forceps



Micro Stan RF™ Forceps

IEC-MJ21

Reusable

Suggested for Blepharoplasty and Facelift procedures


3.5" Length

Includes Quick-Connect™ Cable and Video CD

For use with Surgitron® 4.0 Dual RF™ Device

MJ21

For EMC Device

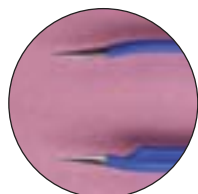

Magnified view



Mini Mono Stan RF™ Forceps

IEC-MJ11/M

Reusable

Monopolar Forceps with Cutting, Coagulation, Irrigation and Suction Port

Suggested for Facelifts


4.875" Length

Includes Quick-Connect™ Cable and Video CD


For use with Surgitron® 4.0 Dual RF™ Device

MJ11/M

For EMC Device

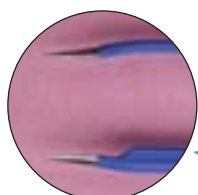

Magnified view

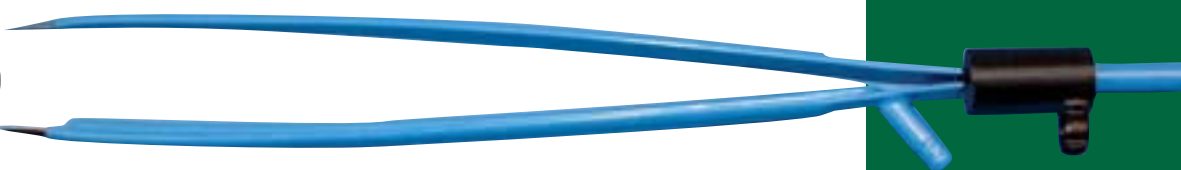

Mono Stan RF™ Forceps

IEC-MJ11

Reusable

Monopolar Forceps with Cutting, Coagulation, Irrigation and Suction Port

Suggested for Breast Augmentation Procedures

6.125" Length

Includes Quick-Connect™ Cable and Video CD

For use with Surgitron® 4.0 Dual RF™ Device

MJ11

For EMC Device

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

The **ellman®** Micro, Mini and Mono Stan Forceps simultaneously dissect and coagulate tissues using the 4.0 MHz radiofrequency.



Facelift

The insides of tips on the Mini Mono and Mono Stan Forceps are serrated.

MULTI-SPECIALTY ACCESSORIES

Dermatology

Cosmetic Treatment Package

For treating telangiectasia and epilation of hair follicles

H135A For EMC Device

IEC-H135A For Surgitron® 4.0 Dual RF™ and Radiolase® II Devices Package Includes:

One MicroInsulated Needle Handpiece

One package of 50 MicroInsulated Sterile Needles and

One package of 50 Patient Education Brochures

XH136 MicroInsulated Needle Handpiece™
For EMC Device (2 pcs)

IEC-XH136 MicroInsulated Needle Handpiece™
For Surgitron® 4.0 Dual™ Devices
(2 pcs)

H137 MicroInsulated Sterile Needles (50 pcs)

H1375 MicroInsulated Sterile Needles – bulk pkg (250 pcs)

LIT/DE1 Patient Education Brochure: FAQ's – Facial Veins (50 pcs)



Magnified view of MicroInsulated Needle

The **ellman**® Cosmetic Treatment Package is an effective, comfortable way to remove small unsightly “spider veins” with minimal recovery time, minimal scarring and immediate improvement in appearance.

With the advent of **ellman**® radiowave sclerotherapy, many small facial veins can now be treated without surgery or injections.¹

Dermatology Electrode Assortment Set

DDERM40

Assorted, Sterile/Single-Use Electrodes, 1/16" shaft (4 electrode styles)



D8D (15 pcs)

B1D (5 pcs)

F1D (15 pcs)

A8D (5 pcs)

Patented Insulated Sclerotherapy™ Needles

D41 27 Gauge, 1" Needle

D42 30 Gauge, 1/2" Needle

D43 Insulated Sclerotherapy Needle Set, 27-Gauge and 30-Gauge
(2 pcs, one of each)



Rhinophyma Electrode Set see page 30

Corneal Shields™ see page 55

Chedly Eyelid Retractor see page 55

For orders and/or pricing information, please call **ellman**® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com



Telangiectasia
Before



Telangiectasia
After

1. Pfenninger, J et al, *Procedures for Primary Care* (2003); second Edition, p221

MULTI-SPECIALTY ACCESSORIES

Non Ablative Skin Tightening Handpieces

FDA 510(k) pending. Skin Tightening product is only available to physicians outside the United States.

Tighter skin without surgery

Unlike prior generations, today's adults are proactive about their health and appearance. They want to look as young as they feel. With increasing numbers, our younger, more active patients are having the early signs of aging such as "crows feet" and frown lines. They are requesting non-invasive, non-surgical procedures that have little, if any, down time.

DHPK1

Non-Ablative Skin Tightening Handpiece Introductory Kit, Reusable.

Includes one each of the 5mm, 10mm, 15mm, and 20mm Dome Handpieces.
One package of Corneal Eye Shields, 25 patient education brochures and instructions for use.
Only used with the Surgitron® Dual Frequency RF™/120 IEC for skin tightening.



DHP05R

Non-Ablative Skin Tightening Handpiece, 5mm limited reusable



DHP10R

Non-Ablative Skin Tightening Handpiece, 10mm limited reusable



DHP15R

Non-Ablative Skin Tightening Handpiece, 15mm limited reusable



DHP20R

Non-Ablative Skin Tightening Handpiece, 20mm limited reusable

(See page 36 for Skin Tightening Patient Education Brochure, Cat. No. LIT/DE2
(See page 55 for Corneal Shields, Cat. No. H781

**For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com**



Close up of tapered-head electrode

MULTI-SPECIALTY ACCESSORIES

Ophthalmology/Oculoplastics

Goisis RF Elevator ACe-Tip™

GE1

Reusable
(1 pc)



Patented RF Wedge Resection Electrode Set™

EE255

Reusable
(set of three)

Patented RF Wedge Resection Electrode

TEE237 1/4" tip (2 pcs)



TEE236 3/16" tip (2 pcs)



TEE235 1/8" tip (2 pcs)



Beigi Minimally Invasive Entropion Electrode

TEE603

Sterile/Single-Use
(2 pcs)



Ophthalmic Micro-Dissection Electrodes™

TA10/C1


Reusable. 2" Electrodes
1.5mm of active .004 Wire (2 pcs)

TA10/C2

Reusable. 2" Electrodes
3.0mm of active .004 Wire (2 pcs)

TA10/C3

Reusable. 2" Electrodes
5.0mm of active .004 Wire (2 pcs)

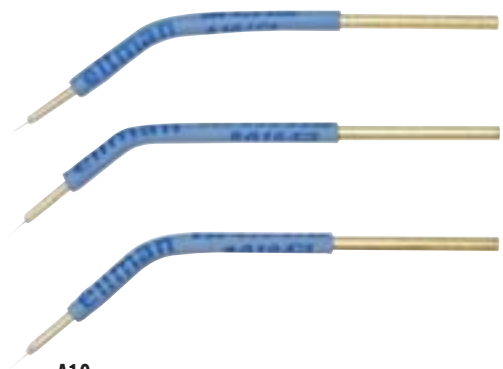

A10

Reusable
Set of two (1.5mm and 3.0mm, .004 Wire - one of each)

Note: Add/OP after the Catalog Number to order 1/2" shorter,
Ophthalmic Micro-Dissection Electrodes™



Goisis RF Elevator ACe-Tip™

Product Benefits:

- Aesthetic midface lift easily done under local anaesthesia.
- Easier and faster orbital decompression in Basedow disease.
- Easier and faster treatment of orbital and facial traumas.
- Versatility – evisceration of the eye with a single instrument.
- Conjunctival peritomy
- Cornea excision
- Globe content removal and full-thickness sclerotomies.

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

MULTI-SPECIALTY ACCESSORIES

Ophthalmology/Oculoplastics

Javate-Pamintuan Dacryoplasty (JPD) Electrode

TEE606

Reusable

(2 pcs)



Corneal Shields™

H781

Autoclavable

Assorted Package
one pair each size,
small, medium, large and
1 suction handle



H75A Small - Red (1 handle, 2 pair)

H76A Medium - Green (1 handle, 2 pair)

H77A Large - Blue (1 handle, 2 pair)

H78A Suction Handle (2 pcs)

H78B Suction Handle (25 pcs)



Chedly Eyelid Retractor

H21A

Reusable

(1 pc)



Baker Brow Probe

EE279

Reusable

(2 pcs)



Magnified view



The **ellman®** Corneal Shield™ tightens the skin to allow for more precise incisions. It is an excellent aid when performing periorbital surgery and aids in the protection from cautery, desiccation, irritation or trauma, during and after surgery of eyelid and orbital regions.

The Chedly Eyelid Retractor allows access for oculoplastic procedures i.e. upper and lower blepharoplasty, ptosis, trichiasis and minor eyelid procedures.

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

MULTI-SPECIALTY ACCESSORIES

Family Practice/Podiatry



ellman® Digit Tourniquets
H97 ellman® Cots – Small (12 pcs)



H98 ellman® Cots – Medium (12 pcs)



H99 ellman® Cots – Large (12 pcs)

ECOTS- 4 each small, medium and large ellman® Cots.



Patented Insulated Matrixectomy Electrodes™*

H9 Reusable
Set of 3 electrodes
1/16" or 3/32" shaft
Matrixectomy Electrode Set
Includes one each of: H10D, H10E, H10A



H10D 2mm Matrix Electrodes
Reusable (2 pc)



H10E 4mm Matrix Spicule Electrodes
Reusable (2 pc)



H10A Matrix Spicule (Ball) Electrode
Reusable (1 pc)



Patented Standard Podiatry Electrode Set*

S11B Reusable
1/16" or 3/32" shaft
Includes one each of: D8, A3, B4, B5,
F1, A8
(6 pcs)



**Electrodes are displayed in actual size.*



Ingrown Toenail

The ellman® radiosurgical nail matrixectomy technique for Ingrown Toenail allows the matrix to be destroyed electrically with a teflon coated electrode, which protects the overlying eponychium.

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

MULTI-SPECIALTY ACCESSORIES

Hemostasis Electrode

Product Benefits

- Quick coagulation can be achieved with our new tapered-head electrode designed for larger area hemostasis
- Available as a 3-Button Fingerswitch Activated Handpiece

3-Button Fingerswitch Activated Handpiece with Integrated Electrode



D110FSHP Hemostasis Electrode 3-Button Fingerswitch Handpiece, 10mm

For orders and/or pricing information, please call ellman® at: (800) 835-5355
(516) 594-3333 • fax (516) 569-0054 • www.ellman.com

INDEX

Accessories	36	Assorted Needle Electrode Set	33	IEC Neutral Plate	9
Educational Accessories	36	Assorted Oval Loop Electrode Set	34	IEC Disposable Neutral Plate-Small Size	9
<i>(Educational Materials)</i>		Assorted Round Loop Electrode Set	33	IEC Disposable Neutral Plate-Standard Size	9
Compact Discs		Assorted Triangle Loops Electrode Set	34	Non Ablative Skin Tightening Handpieces	53
Oculoplastic Techniques		Rodin® Microsurgery ACe-Tip™ Electrode Set	33	Ophthalmology/Oculoplastics	54
Patient Educational Brochures		Standard Electrode Set	33	Baker Brow Probe	55
Textbooks		Energy Sources	5	Beigi Minimally Invasive Entropion Electrode	54
Device/Accessories	37	Surgitron® Dual RF™	5	Cheddy Eyelid Retractor	55
<i>(Smoke Evacuation System)</i>		Surgitron® FFPF EMC	8	Corneal Shields™	55
SugiCart™		Surgitron® Radiolase® I	6	Goisis RF Elevator ACe-Tip™	54
Surg-e-Vac™		Surgitron® Radiolase® II	7	Javate-Pamintuan Dacryoplasty (JPD) Electrode	55
Vac Assist Stand		Family Practice/Podiatry	56	Ophthalmic Micro-Dissection Electrodes™	54
Accessories	38	ellman® Digit Tourniquets	56	Patented RF Wedge Resection Electrode	54
<i>(Smoke Evacuation System)</i>		Patented Insulated Matrixectomy Electrodes™	56	Patented RF Wedge Resection Electrode Set™	54
Surg-e-Vac™ Accessories		Patented Standard Podiatry Electrode Set	56	Otorhinolaryngology	40
For both Surg-e-Vac™ and Vapor-Vac II™ Accessories		Footswitches	10	3" Length Handpiece with Patented Locking Chuck	40
Accessories	39	Dual Frequency Footswitch and Cable Pneumatic Footswitch	10	Adenoid RF Electrode™	43
<i>(Multi-Specialty Accessories)</i>		Radiolase I Footswitch and Cable	10	Bayonet RF Turbinate Volume Reduction Electrode™	42
Medi-Frig™		Radiolase II Monopolar/Bipolar Footswitch	10	Cleaning Brush for Suction Probe	44
Wipe-Its™		General Surgery/Proctology	47	Clear-Vu® Bayonet Electrodes	42
ACe-Tip™ Low-Stick Forceps	20-21	Patented Gupta RF Gun (Pravin Gupta Procto-Gun)	47	Cryptolysis Electrode	46
Adapters For Surgitron® EMC	13	Patented RF Laparoscopic Electrodes™	47	HLTZ Tonsillotomy Electrode™	43
Bipolar Converter Steri Caps™	14	Patented RF Laparoscopic Handpiece	47	Micro-Larynx™ RF Probes	40
EMC Steri Caps™	14	Gynecology	48	OTO-RAM™ Radiofrequency Assisted Myringotomy Set	41
Patented 1-Button Fingerswitch Handpiece with Locking Chuck	13	Conization Electrode™	48	Patented Micro Fibre® Electrode	40
Patented Bipolar Converter	14	External Genital Lesion RF Electrodes™	48	Patented OTO-RAM™ Probes	41
Universal Fingerswitch Adapter	13	LLETZ Bayonet Electrodes	49	Patented RAUP Electrode™	43
Wall Bracket	14	LLETZ Electrodes Assortment Set	48	Patented RF Angled Suction Coagulator™	45
Bayonet Bipolar Forceps	19	LLETZ Individual Electrodes	48	Patented RF Bipolar Tonsil Reduction Probe™ and Quick-Connect™ Cable	44
Bipolar Forceps	17-18	Monssel's Solution	49	Patented RF Myringotomy Electrode™	41
Cables	11	Needle Extender/Measured Dose Syringe	49	Patented RF Suction Adenoid Probes	44
Bipolar Cable	11	Standard Electrode Set	33	Patented RF Suction Coagulator™	45
Lap/Endo Monopolar Cable	11	Vaginal Lateral Wall Retractor	49	Patented Suction Coagulation Kit	45
Snare Cable	11	Vaginal Speculums	49	Patented Suction Tonsillectomy Probe™ and Quick-Connect™ Cable	44
Quick-Connect™ Cable	11	Handpieces	15	Patented Tongue Base Reduction Electrode™	43
Dermatology	52	Blade Handpiece with Locking Chuck	16	Polyp Snare™ System for Nasal Polyps	46
Cosmetic Treatment Package	52	Foot-Controlled Handpiece with Locking Chuck	16	Quick-Connect™ Cable	40
Dermatology Electrode Assortment Set	52	microlnsulated Needle Handpiece with Locking Chuck	16	RF Assisted Tonsillotomy (RAT) Ball Electrodes™	46
Patented Insulated Sclerotherapy™ Needles	52	One-Button Fingerswitch Handpiece	16	Sub-Mucosal VTR (UPP) Electrode	43
Electrodes	23	Patented 1-Button Fingerswitch Handpiece with Locking Chuck	15	Plastic Surgery	50
Electrodes Accessories	35	Patented 3-Button Fingerswitch Blade Handpiece with Locking Chuck	15	Javate/Khan Coated Flat Scalpel	50
Sterile/Single-Use	24	Patented 3-Button Fingerswitch Handpiece with Locking Chuck	15	Javate/Khan Coated Spade	50
Empire microlnscision™ Needle Electrodes	25-26	Rodin® Microsurgery Handpiece	16	Suction Set	50
Medical Electrode Set	24	Hemostasis Electrode	57	Micro Stan RF™ Forceps	51
Reusable Electrodes	27	3-Button Fingerswitch Activated	57	Mini Mono Stan RF™ Forceps	51
ACe-Tip™		Neutral Plates	9	Mono Stan RF™ Forceps	51
Advanced Composition Alloy	32	EMC Neutral Plate	9	RF SurgiPen™	22
Ball Electrodes	29			Surg-e-Doc™	12
Diamond Loop Electrodes	28				
Empire® MicroTip Electrodes	30				
Insulated Needle Electrodes™	30				
Javate Endoscopic DCR Electrodes	31				
Needle Electrodes	27				
Oval Loop Electrodes	29				
Radiosurgical Mucotome Electrodes	31				
Rhinophyma Electrode Set	30				
Round Loop Electrodes	27				
Scalpel Electrode	30				
Triangle Loop Electrodes	28				
Reusable Electrodes Sets	33				
Assorted Ball Electrode Set	34				
Assorted Diamond Loop Electrode Set	34				

ellman®

Experts in Precision Surgery

Established in 1959



ellman® International, Inc.

3333 Royal Avenue, Oceanside, N.Y. 11572-3625 U.S.A.
(800) 835-5355 • (516) 594-3333 • Fax: (516) 569-0054
www.ellman.com

CC09003SVF

Copyright 2009 ellman® International, Inc. All Rights Reserved.

Note: Do not copy or distribute without written authorization from Ellman International Inc.

